### **Features**

Supports the SensL ArrayC-30035-144P-PCB 12x12 array of 3mm SiPMs

Horizontal signal connector located on the back, array located on the front

4-side tileable installation

Four encoded position signals for event centroid calculations: X+, X-, Y+, Y-

DC-coupled signal path

Low power consumption

Patented diode-coupled charge division readout, superior to traditional resistive readout

Improved spatial uniformity

Faster rise time


Reduced image noise

Precision temperature sensor

Four mounting holes for #2 hardware

SensL's fast output signals are not connected



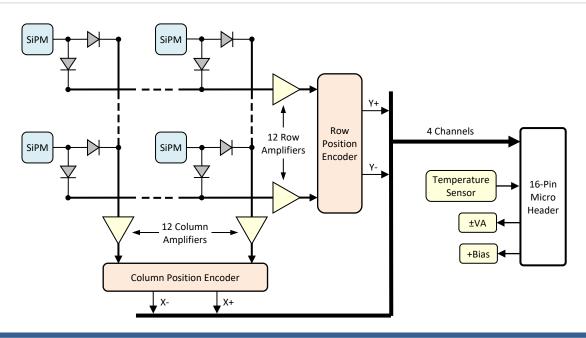

# **Specifications**

### **Position Signal Outputs**

Encoding Charge division multiplexed to

4 output channels: X+, X-, Y+, Y-

Gain  $750\Omega$  transimpedance gain

Output voltage  $0 \rightarrow -1V$  into  $100\Omega$  load

Output impedance  $100\Omega$ 

Output current 50mA maximum

**Temperature Sensor** 

Output voltage 500mV + 10mV per °C

Output current 10mAOutput impedance  $100\Omega$ Accuracy  $\pm 0.5^{\circ}C$ 

Bias Voltage +29V typical (refer to SiPM data)

Voltage clamp 47V Zener diode

500mW maximum

Amplifier Voltage ( $\pm$ VA)  $\pm 2.8V \rightarrow \pm 5.5V$  DC maximum

Current ±50mA typical

(Iq, no signal, no load)


Signal Connector Horizontal 16-pin 2-row header

0.050" pin pitch


Mating assembly Samtec FFSD-08-D-XX.XX-01-N

(XX.XX = length in inches)

### **Channel Map**

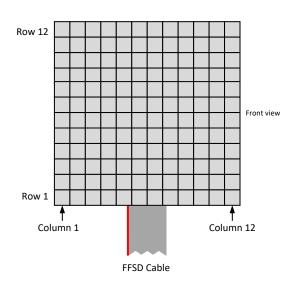

### **Signal Connector**

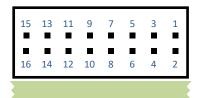

Side View

Pin	Function	Pin	Function
1	Temperature	2	Ground
3	X-	4	Ground
5	X+	6	Ground
7	-VA	8	Ground
9	+VA	10	Ground
11	Υ-	12	Ground
13	Y+	14	Ground
15	+Bias	16	Ground

### **4-Channel Position Encoder**

### **Row and Column Encoder Weights**

Row# or Col# (for X- or Y-)	Row# or Col# (for X+ or Y+)	Fraction (ideal)	Fraction (actual)	% Error	Notes
1	12	0.0833	0.0833	0.00 %	
2	11	0.1667	0.1650	-1.02 %	
3	10	0.2500	0.2483	-0.68 %	
4	9	0.3333	0.3311	-0.66 %	
5	8	0.4167	0.4167	0.00 %	Sum of X- and X+ fractions
6	7	0.5000	0.5000	0.00 %	or Y- and Y+ fractions
7	6	0.5833	0.5882	0.84 %	= 1.0833
8	5	0.6667	0.6637	-0.45 %	Independent of signal position
9	4	0.7500	0.7500	0.00 %	
10	3	0.8333	0.8333	0.00 %	
11	2	0.9167	0.9091	-0.83 %	
12	1	1.0000	1.0000	0.00 %	

**Note:** Errors exclude component tolerances

### **Output Signals**

X- = (SiPM signal) \* (encoder gain) \* (X- fraction)

X+ = (SiPM signal) \* (encoder gain) \* (X+ fraction)

Y- = (SiPM signal) \* (encoder gain) \* (Y- fraction)

Y+ = (SiPM signal) \* (encoder gain) \* (Y+ fraction)

Typical event position calculation:

X column = (X + - X -) / (X + + X -)

Y row = (Y+ - Y-) / (Y+ + Y-)

### **Example**

SiPM signal at column 4, row 3 (excluding encoder gain)

X- = (Column 4 signal) \* 0.3311

**X+** = (Column 4 signal) \* 0.7500


**Y-** = (Row 3 signal) \* 0.2483


**Y+** = (Row 3 signal) \* 0.8333


### Datasheet Rev. A-1811




# **Typical Signals**

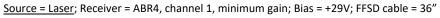

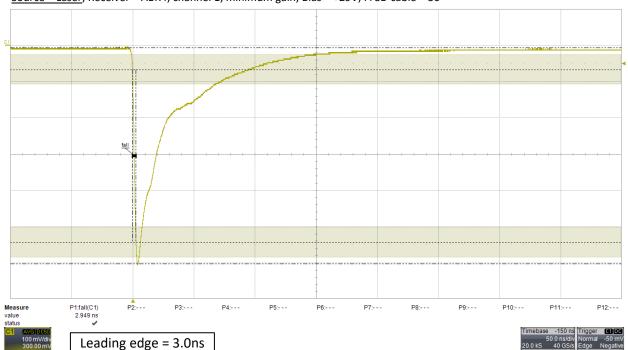

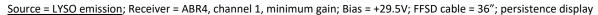


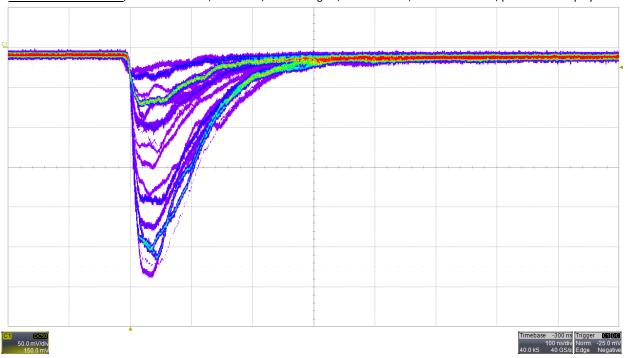



### Mechanical


### Front View Back View

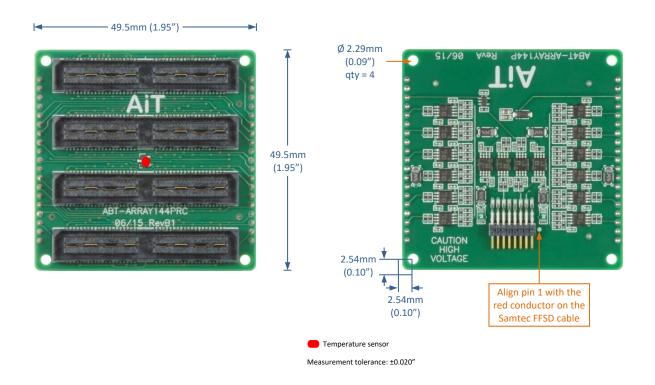


### Side View, Base Attached to Array

# Base 18.1mm (0.725") 3.56mm (0.14")

### Side View, Base Only

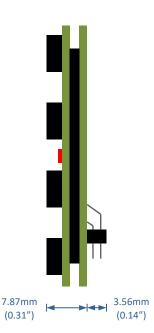

### 4-Channel Active Base Readout Kit

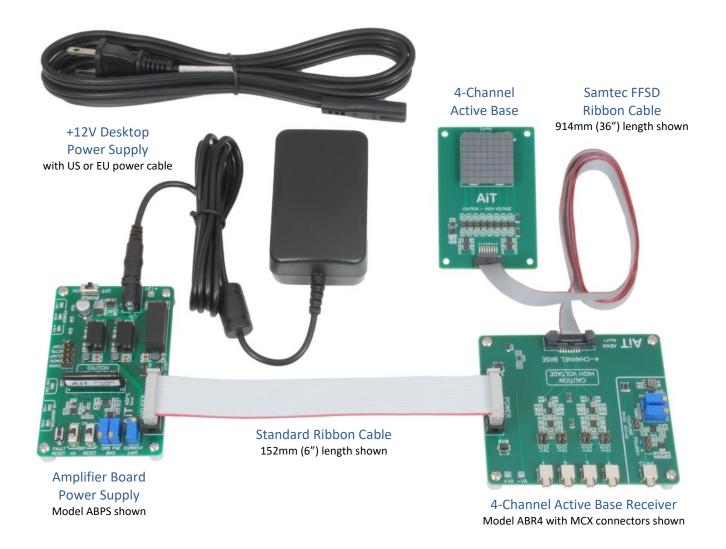

### **Components**

Each component is available separately. Refer to each datasheet for details.

The Active Base includes a 914mm (36") Samtec FFSD micro-pitch ribbon cable.

The Amplifier Board Power Supply includes a 12V desktop power supply and a HV80 bias voltage power supply.

The 4-channel Active Base Receiver includes a 152mm (6") power supply ribbon cable and a breakout board to connect any external power supply.

# **Safety Information**



# WARNING – High Voltage

- High voltage may be present during operation
- High voltage stored on capacitors may be present after power is removed
- Improper handling may result in personnel injury or equipment damage

This high-voltage device must be used only by personnel trained and qualified in safe handling, installation, and operation of high-voltage equipment.



# **CAUTION – Electrostatic Discharge (ESD) Sensitivity**

The circuit board can be damaged by electrostatic discharge. Observe precautions for handling electrostatic sensitive devices. Handle only at static-safe workstations.

## **High-Gain Photodetectors**

High-gain photodetectors such as silicon photomultipliers may conduct damaging currents if exposed to high optical signal levels while the bias voltage is applied, or if the bias voltage exceeds the recommended operating range. These devices must be operated only in low-light conditions, and only within the manufacturer's recommended bias voltage range.

### **Handling and Disassembly**

This product may be provided with a protective enclosure. Disassembled enclosure components and circuit boards may contain sharp edges. Take appropriate safety precautions while assembling or disassembling the enclosure and handling disassembled components.

### **Indoor Use Only**

Do not operate this product in a wet or damp environment. Do not operate in an explosive atmosphere.

Use of this product, and AiT Instruments' liability related to use of this product, is further governed by AiT Instruments' standard terms and conditions of sale, which were provided upon purchase of this product.