

High-Speed Balanced Photoreceiver




Features

- **Bandwidth DC to 500 MHz**
- **Common-Mode Rejection Ratio (CMRR) 40 dB typ.**
- **SI-PIN detectors, 0.4 mm active diameter**
- **Spectral range 320 - 1000 nm**
- **Very low NEP, down to 12 pW/√Hz**
- **Transimpedance gain switchable 5×10^3 V/A, 10×10^3 V/A**
- **High dynamic input range up to 2×10 mW balanced optical power**
- **Fast monitor outputs with 10 MHz bandwidth and 1×10^3 V/A gain**
- **Switchable low pass filter for minimizing wideband noise**
- **Free-space input 1.035"-40 threaded, alternatively 25 mm diameter unthreaded**
- **Easily convertible to fiber optic input (FC and FSMA) with optionally available screw-on adapters**
- **UNC 8-32 and M4 tapped holes for mounting on standard posts with metric and imperial thread**

Applications

- **Spectroscopy**
- **Heterodyne detection**
- **Optical coherence tomography (OCT)**
- **Optical delay measurement**
- **Differential optical front-end for oscilloscopes, spectrum analyzers, A/D converters and RF lock-in amplifiers**

Block Diagram

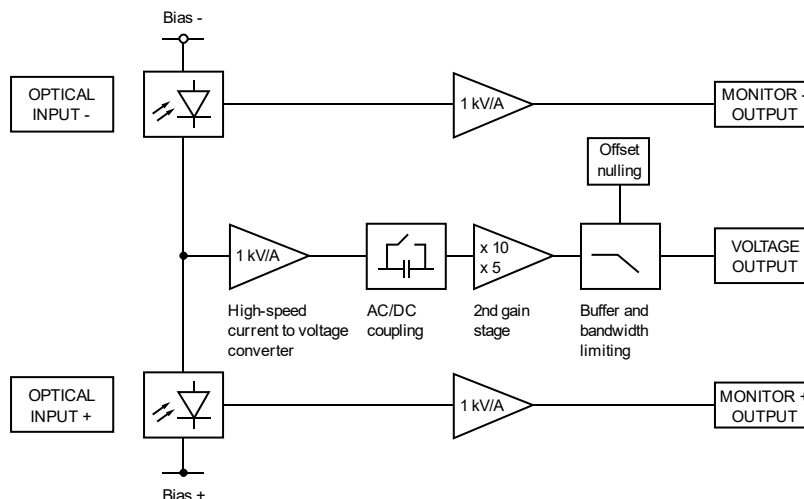

High-Speed Balanced Photoreceiver

Available Input Versions

HBPR-500M-10K-SI-FST



Picture shows two 1.035"-40 threaded flanges with internally threaded coupler rings mounted (outer diameter 30 mm)

1.035"-40 threaded flange for free space applications, compatible with many optical standard accessories.

Optional: Fiber adapters PRA-FC, PRA-FCA, PRA-FSMA



HBPR-500M-10K-SI-FS



25 mm dia. unthreaded flange for free space applications compatible with many optical standard accessories.

Related Models

Various free space or fiber coupled HBPR models, with bandwidth up to 500 MHz, in the spectral range from 320 nm to 1700 nm are available.

Example: FC input



fix/permanent FC fiber connector for high coupling efficiency, excellent conversion gain accuracy and common mode rejection ratio (CMRR).

See further information and separate datasheets on www.femto.de

Available Accessories

PRA-FC
PRA-FCA
PRA-FSMA



fiber-adapter with external 1.035"-40 thread (suitable for FST models only)

PS-15




power supply,
input: 100 - 240 VAC,
output: ± 15 VDC, +400/-250 mA

High-Speed Balanced Photoreceiver


Specifications	Test conditions	$V_S = \pm 15\text{ V}$, $T_A = 25\text{ }^\circ\text{C}$, signal output terminated with $50\text{ }\Omega$, Monitor outputs terminated with $1\text{ M}\Omega$
Gain	Transimpedance gain	$5 \times 10^3\text{ V/A}$ (2 nd gain x5), $10 \times 10^3\text{ V/A}$ (2 nd gain x10) switchable (@ $50\text{ }\Omega$ load)
	Gain accuracy	$\pm 1\%$ electrical
	Conversion gain	$2.55 \times 10^3\text{ V/W}$ typ. (@ 2 nd gain x5, 760 nm) $5.1 \times 10^3\text{ V/W}$ typ. (@ 2 nd gain x10, 760 nm)
	Common mode rejection ratio (CMRR)	50 dB typ. ($f \leq 100\text{ MHz}$) 40 dB typ. ($f \leq 500\text{ MHz}$)
Frequency Response	Lower cut-off frequency	DC / 10 Hz, switchable
	Upper cut-off frequency	500 MHz (@ 2 nd gain x5); 460 MHz (@ 2 nd gain x10), switchable to 20 MHz
Time Response	Rise/fall time (10 % - 90 %)	0.85 ns (@ 2 nd gain x5); 0.95 ns (@ 2 nd gain x10) 17.5 ns (low pass filter 20 MHz)
Input	Noise equivalent power (NEP)	minimum $12\text{ pW}/\sqrt{\text{Hz}}$ (@ 760 nm) $13\text{ pW}/\sqrt{\text{Hz}}$ (@ 760 nm, 20 MHz) $29\text{ pW}/\sqrt{\text{Hz}}$ (@ 760 nm, 200 MHz) $60\text{ pW}/\sqrt{\text{Hz}}$ (@ 760 nm, 500 MHz)
	Maximum differential CW power for linear amplification	$400\text{ }\mu\text{W}$ (@ 2 nd gain x5, DC-coupled, 760 nm) $200\text{ }\mu\text{W}$ (@ 2 nd gain x10, DC-coupled, 760 nm) 2.5 mW (@ AC-coupled, 760 nm)
	Max. optical CW balanced power (common mode power)	10 mW (on each photodiode, @ 760 nm)
	Monitor optical saturation power (limited by Maximum Rating)	12 mW (@ 760 nm)
Detector	Detector	SI-PIN photodiode
	Active area	$\varnothing 400\text{ }\mu\text{m}$
	Spectral range	320 - 1000 nm
	Sensitivity	0.51 A/W typ. (@ 760 nm)
Signal Output	Output voltage range	$\pm 1.0\text{ V}$ (@ $50\text{ }\Omega$ load) for linear operation and low harmonic distortion
	Max. output voltage	$\pm 2.0\text{ V}$ (@ $50\text{ }\Omega$ load)
	Offset voltage compensation	$\pm 100\text{ mV}$ typ., adjustable by offset potentiometer
	Output impedance	$50\text{ }\Omega$ (terminate with $50\text{ }\Omega$ load)
	Slew rate	2800 V/ μs
	Max. output current	70 mA
	Output return loss S22	-30 dB @ $< 100\text{ MHz}$ -20 dB @ $< 800\text{ MHz}$
	Output noise	$2.5\text{ mV}_{\text{RMS}}$ (16 mV _{PP}) (@ 2 nd gain x5) $4.1\text{ mV}_{\text{RMS}}$ (27 mV _{PP}) (@ 2 nd gain x10) $0.25\text{ mV}_{\text{RMS}}$ (1.7 mV _{PP}) typ. (@ 2 nd gain x5, BW: 20 MHz) $0.4\text{ mV}_{\text{RMS}}$ (2.5 mV _{PP}) typ. (@ 2 nd gain x10, BW: 20 MHz) (@ $50\text{ }\Omega$ load, no signal on detectors, measurement bandwidth 2 GHz)

High-Speed Balanced Photoreceiver

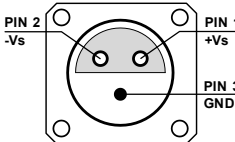Specifications (Continued)																			
Monitor Outputs	<table><tr><td>Monitor output gain</td><td>1 x 10³ V/A (@ ≥ 100 kΩ load)</td></tr><tr><td>Monitor output voltage range</td><td>0 ... +10 V (@ ≥ 100 kΩ load)</td></tr><tr><td>Monitor output impedance</td><td>50 Ω (terminate with ≥ 100 kΩ load)</td></tr><tr><td>Monitor output max. output current</td><td>30 mA typ.</td></tr><tr><td>Monitor output bandwidth</td><td>DC ... 10 MHz</td></tr><tr><td>Monitor output noise</td><td>0.6 mV_{RMS} (4 mV_{PP}) (@ 100 kΩ load, no signal on detectors, measurement bandwidth 200 MHz)</td></tr></table>	Monitor output gain	1 x 10 ³ V/A (@ ≥ 100 kΩ load)	Monitor output voltage range	0 ... +10 V (@ ≥ 100 kΩ load)	Monitor output impedance	50 Ω (terminate with ≥ 100 kΩ load)	Monitor output max. output current	30 mA typ.	Monitor output bandwidth	DC ... 10 MHz	Monitor output noise	0.6 mV _{RMS} (4 mV _{PP}) (@ 100 kΩ load, no signal on detectors, measurement bandwidth 200 MHz)						
Monitor output gain	1 x 10 ³ V/A (@ ≥ 100 kΩ load)																		
Monitor output voltage range	0 ... +10 V (@ ≥ 100 kΩ load)																		
Monitor output impedance	50 Ω (terminate with ≥ 100 kΩ load)																		
Monitor output max. output current	30 mA typ.																		
Monitor output bandwidth	DC ... 10 MHz																		
Monitor output noise	0.6 mV _{RMS} (4 mV _{PP}) (@ 100 kΩ load, no signal on detectors, measurement bandwidth 200 MHz)																		
Input Flange	<table><tr><td>Material</td><td>1.4305 stainless steel, nickel-plated (FST flange) AlMg4.5Mn, nickel-plated (FS flange)</td></tr></table>	Material	1.4305 stainless steel, nickel-plated (FST flange) AlMg4.5Mn, nickel-plated (FS flange)																
Material	1.4305 stainless steel, nickel-plated (FST flange) AlMg4.5Mn, nickel-plated (FS flange)																		
Coupler Ring (FST version only)	<table><tr><td>Material</td><td>1.4305 stainless steel, glass bead blasted</td></tr></table>	Material	1.4305 stainless steel, glass bead blasted																
Material	1.4305 stainless steel, glass bead blasted																		
Power Supply	<table><tr><td>Supply voltage</td><td>±15 V (±14.5 V ... ±16.5 V)</td></tr><tr><td>Supply current</td><td>−90 / +120 mA (depends on operating conditions, recommended power supply capability min. ±200 mA)</td></tr></table>	Supply voltage	±15 V (±14.5 V ... ±16.5 V)	Supply current	−90 / +120 mA (depends on operating conditions, recommended power supply capability min. ±200 mA)														
Supply voltage	±15 V (±14.5 V ... ±16.5 V)																		
Supply current	−90 / +120 mA (depends on operating conditions, recommended power supply capability min. ±200 mA)																		
Case	<table><tr><td>Weight</td><td>350 g (0.77 lbs)</td></tr><tr><td>Material</td><td>AlMg3Mn, nickel-plated</td></tr></table>	Weight	350 g (0.77 lbs)	Material	AlMg3Mn, nickel-plated														
Weight	350 g (0.77 lbs)																		
Material	AlMg3Mn, nickel-plated																		
Temperature Range	<table><tr><td>Storage temperature</td><td>−40 ... +85 °C</td></tr><tr><td>Operating temperature</td><td>0 ... +60 °C</td></tr></table>	Storage temperature	−40 ... +85 °C	Operating temperature	0 ... +60 °C														
Storage temperature	−40 ... +85 °C																		
Operating temperature	0 ... +60 °C																		
Absolute Maximum Ratings	<table><tr><td>Max. CW power (averaged)</td><td>12 mW (on each photodiode)</td></tr><tr><td>Power supply voltage</td><td>±20 V</td></tr></table>	Max. CW power (averaged)	12 mW (on each photodiode)	Power supply voltage	±20 V														
Max. CW power (averaged)	12 mW (on each photodiode)																		
Power supply voltage	±20 V																		
Connectors	<table><tr><td>Input</td><td>FS version</td><td>25 mm dia. unthreaded flange for free space applications</td></tr><tr><td></td><td>FST version</td><td>1.035"-40 threaded flange for free space applications and for use with various types of optical standard accessories</td></tr><tr><td>Output</td><td>SMA jack (female)</td><td></td></tr><tr><td>Power supply</td><td>Lemo® series 1S, 3-pin fixed socket (mating plug type: FFA.1S.303.CLAC52)</td><td></td></tr></table> <table><tr><td>Pin 1:</td><td>+15 V</td></tr><tr><td>Pin 2:</td><td>−15 V</td></tr><tr><td>Pin 3:</td><td>GND</td></tr></table>	Input	FS version	25 mm dia. unthreaded flange for free space applications		FST version	1.035"-40 threaded flange for free space applications and for use with various types of optical standard accessories	Output	SMA jack (female)		Power supply	Lemo® series 1S, 3-pin fixed socket (mating plug type: FFA.1S.303.CLAC52)		Pin 1:	+15 V	Pin 2:	−15 V	Pin 3:	GND
Input	FS version	25 mm dia. unthreaded flange for free space applications																	
	FST version	1.035"-40 threaded flange for free space applications and for use with various types of optical standard accessories																	
Output	SMA jack (female)																		
Power supply	Lemo® series 1S, 3-pin fixed socket (mating plug type: FFA.1S.303.CLAC52)																		
Pin 1:	+15 V																		
Pin 2:	−15 V																		
Pin 3:	GND																		
Scope of Delivery	HBPR-500M-10K-SI, 2 x threaded coupler ring (FST version only), Lemo® 3-pin connector, 3 x adapter SMA (male) to BNC (female), datasheet																		
Ordering Information	<table><tr><td>HBPR-500M-10K-SI-FS</td><td>25 mm dia. unthreaded flange for free space applications</td></tr><tr><td>HBPR-500M-10K-SI-FST</td><td>1.035"-40 threaded flange for free space applications and for use with various types of optical standard accessories</td></tr></table>	HBPR-500M-10K-SI-FS	25 mm dia. unthreaded flange for free space applications	HBPR-500M-10K-SI-FST	1.035"-40 threaded flange for free space applications and for use with various types of optical standard accessories														
HBPR-500M-10K-SI-FS	25 mm dia. unthreaded flange for free space applications																		
HBPR-500M-10K-SI-FST	1.035"-40 threaded flange for free space applications and for use with various types of optical standard accessories																		


High-Speed Balanced Photoreceiver

Spectral Responsivity

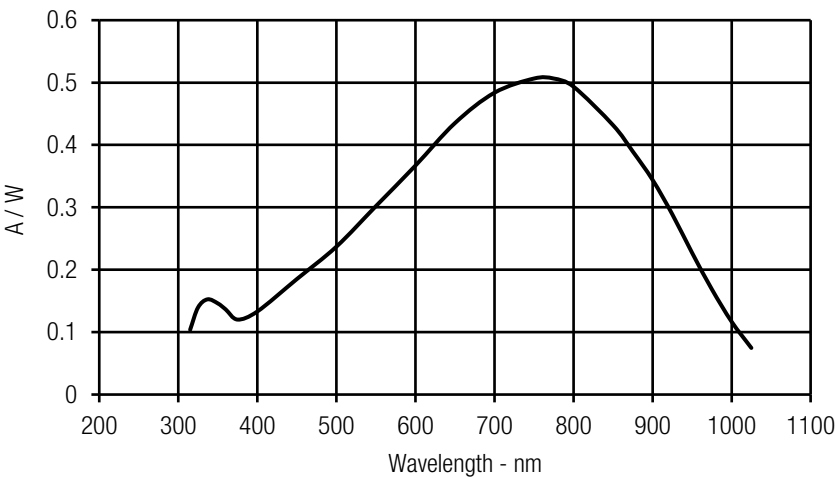


Detector Position

HBPR-500M-10K-SI-FST (1.035"-40 threaded free space input)

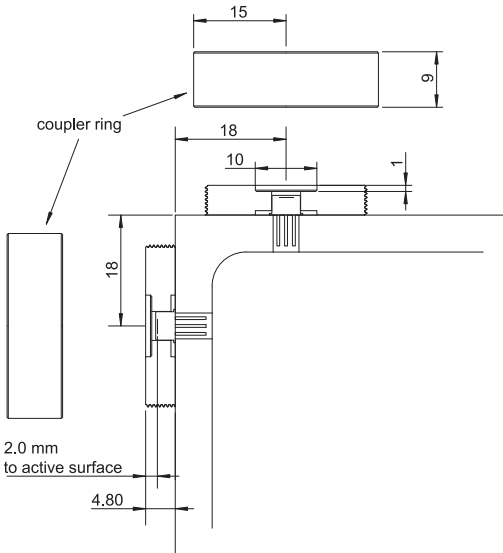


HBPR-500M-10K-SI-FS (25 mm dia. unthreaded free space input)

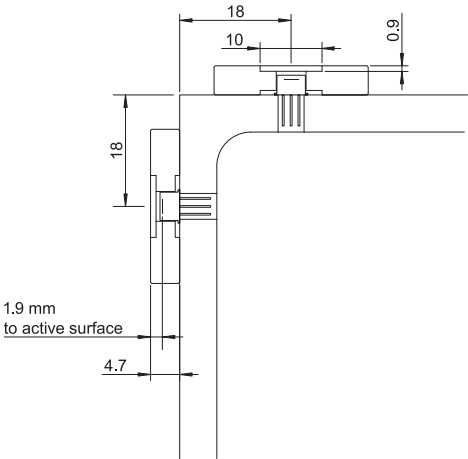

High-Speed Balanced Photoreceiver

Dimensions

Case dimensions for HBPR-500M-10K-SI (FS/FST model):

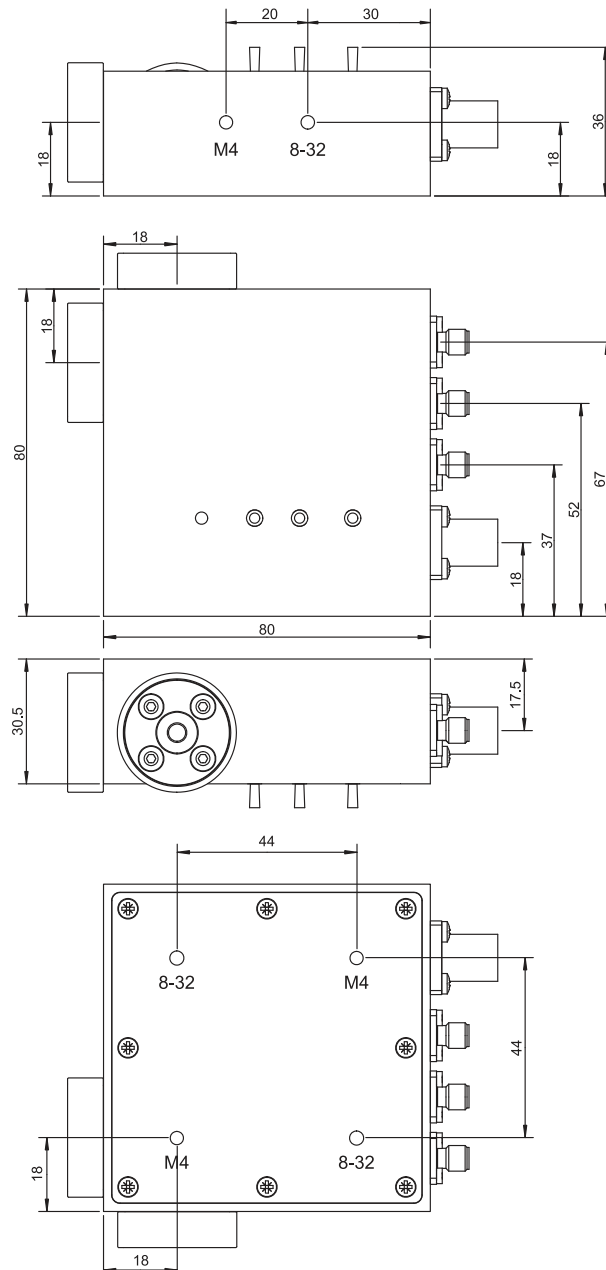

All measures in mm unless otherwise noted.

The bottom plate may be rotated to match the appropriate mounting thread to the optical axis by unscrewing the 8 screws.

FEMTO Messtechnik GmbH
Klosterstr. 64
10179 Berlin · Germany
Phone: +49 30 280 4711-0
Fax: +49 30 280 4711-11
Email: info@femto.de
www.femto.de

Specifications are subject to change without notice. Information provided herein is believed to be accurate and reliable. However, no responsibility is assumed by FEMTO Messtechnik GmbH for its use, nor for any infringement of patents or other rights of third parties which may result from its use. No license is granted by implication or otherwise under any patent or patent rights of FEMTO Messtechnik GmbH. Product names mentioned may also be trademarks used here for identification purposes only.

© by FEMTO Messtechnik GmbH · Printed in Germany